


Vascular Intervention // **Peripheral**
Reinforced Introducer Sheath System

Fortress[®]



Excellent kink and
deformation resistance

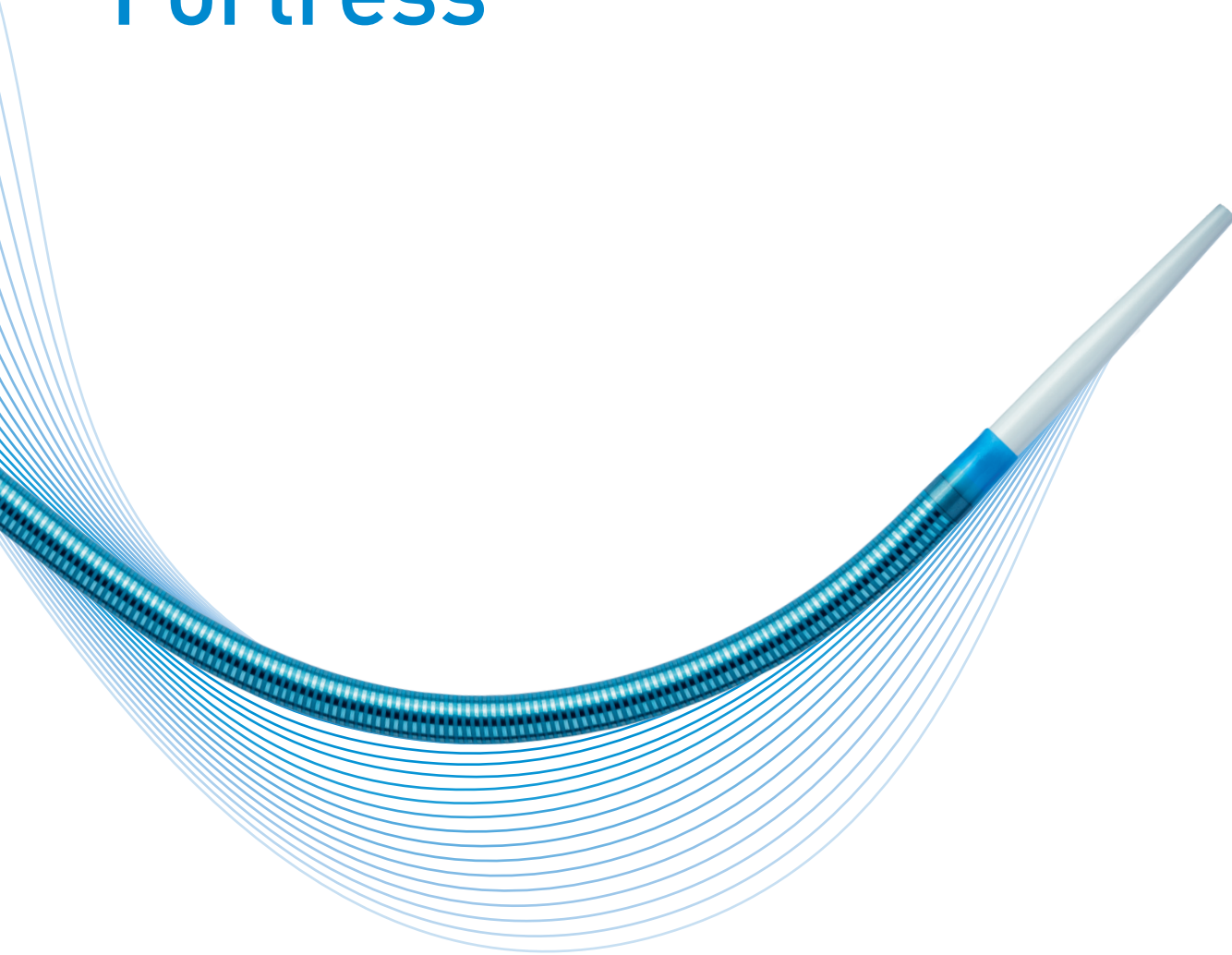

Easy insertion



Superb radiopacity



BIOTRONIK
excellence for life

Fortress

Introducer sheath system with excellent kink resistance, easy insertion and superb radiopacity.


Excellent kink and deformation resistance


Coil reinforced sheath

The polymer embedded stainless steel coil shaft provides flexibility and excellent kink and deformation resistance.

PTFE liner

The non-ovaling characteristic combined with its friction lowering PTFE liner enables devices to easily pass through tight anatomical curves.

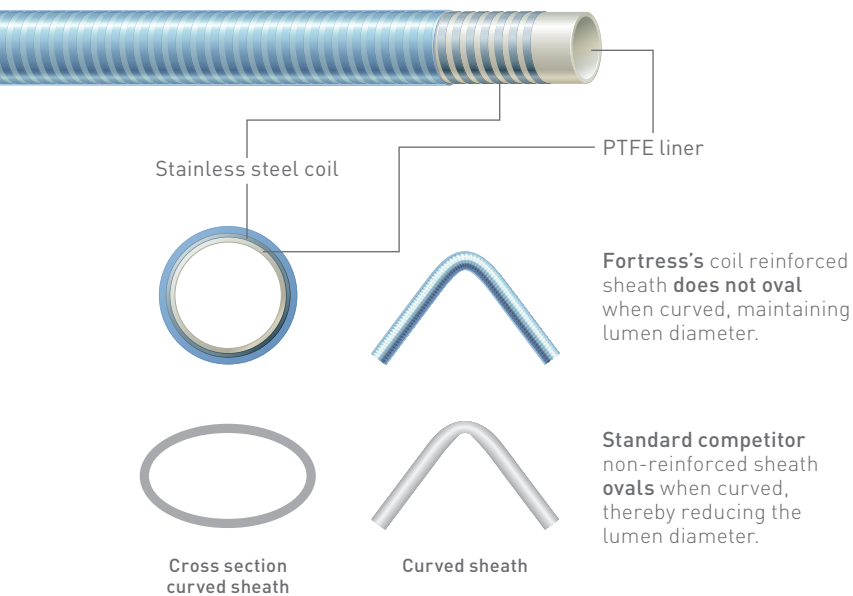

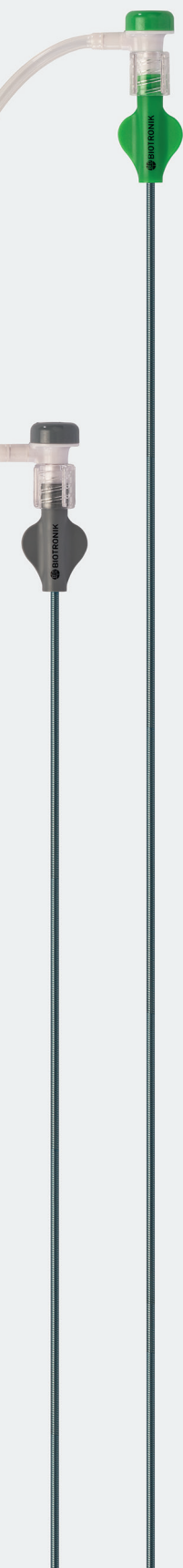

Easy insertion

Smooth tapered sheath and dilator tip

The sheath and dilator have a smooth tapered tip and an atraumatic tip design for ease of insertion.


Hydrophobic coated sheath¹

Best in class insertion performance² due to the hydrophobic sheath coating in combination with the smooth tapered sheath and dilator tip.

Superb radiopacity

Radiopaque sheath marker

The platinum alloy marker band provides superb tip radiopacity under fluoroscopic imaging.³

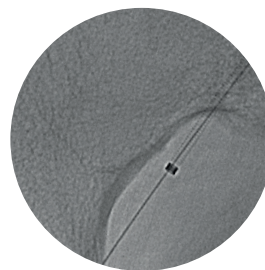

Fortress 4F in the pelvic area

Minimal is optimal

4F

The BIOTRONIK 4F Solutions product portfolio allows for treatment of complex lower limb disease while reducing access site complications.

Fortress

Vascular
Intervention
Peripheral



Indicated to provide access into the femoral, popliteal and infrapopliteal arteries.*

Technical Data

Introducer Sheath

French size	4F, 5F and 6F
Length	45, 65 and 90 cm
Sheath material	Stainless steel coil reinforced polymer, PTFE liner
Sheath tip	Smooth taper with radiopaque marker band located 3-5 mm from the sheath tip
Sheath shape	Pre-curved 45 cm Straight 45, 65 and 90 cm
Sheath coating	Hydrophobic, distal 30 cm (5F and 6F only)

Dilator

Recommended guide wire	0.035"
Tip	Smooth taper with radiopaque filler material
Shape	Pre-curved and straight
Hub	4F: fixed hemostatic valve 5F, 6F: removable hemostatic valve Side-arm with color coded 3-way stop cock

Ordering Information

	Sheath length (cm)	Sheath shape	Introducer description	Catalogue number	Dilator(s) included	Units per box
4F	45	Curved	Crossover	358813	Curved/Straight	5
	45	Straight	Straight-45	358814	Curved/Straight	5
	65	Straight	Straight-65	444486	Straight	5
	90	Straight	Straight-90	444485	Straight	5
5F	45	Curved	Crossover	386590	Straight	5
	45	Straight	Straight-45	386591	Straight	5
	65	Straight	Straight-65	444484	Straight	5
	90	Straight	Straight-90	444483	Straight	5
6F	45	Curved	Crossover	386593	Straight	5
	45	Straight	Straight-45	386594	Straight	5
	65	Straight	Straight-65	444482	Straight	5
	90	Straight	Straight-90	444481	Straight	5

1. 5F and 6F sheath; 2, 3. BIOTRONIK data on file.

*Indication as per IFU.

Fortress is a trademark or registered trademark of the BIOTRONIK Group of Companies.

Manufacturer:

Contract Medical
International GmbH
Lauensteiner Strasse 37
01277 Dresden, Germany
Tel +49 351 213 88 88
Fax +49 351 21388 99
info@contract-medical.com
www.contract-medical.com

Distributor:

BIOTRONIK AG
Ackerstrasse 6
8180 Bülach, Switzerland
Tel +41 (0) 44 8645111
Fax +41 (0) 44 8645005
info.vi@biotronik.com
www.biotronik.com

© 2019 BIOTRONIK AG – All rights reserved.
Specifications are subject to modification, revision and improvement.

 **BIOTRONIK**
excellence for life