


Up to 3.8 x faster
deflation times



Enhanced crossability

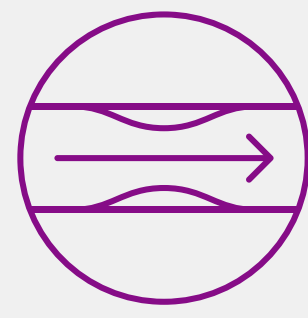

High pushability
and flexibility




Technical data /
ordering info

Vascular Intervention // **Peripheral**
PTA Balloon Catheter/0.014"/OTW

 **BIOTRONIK**
excellence for life

Passeo-14

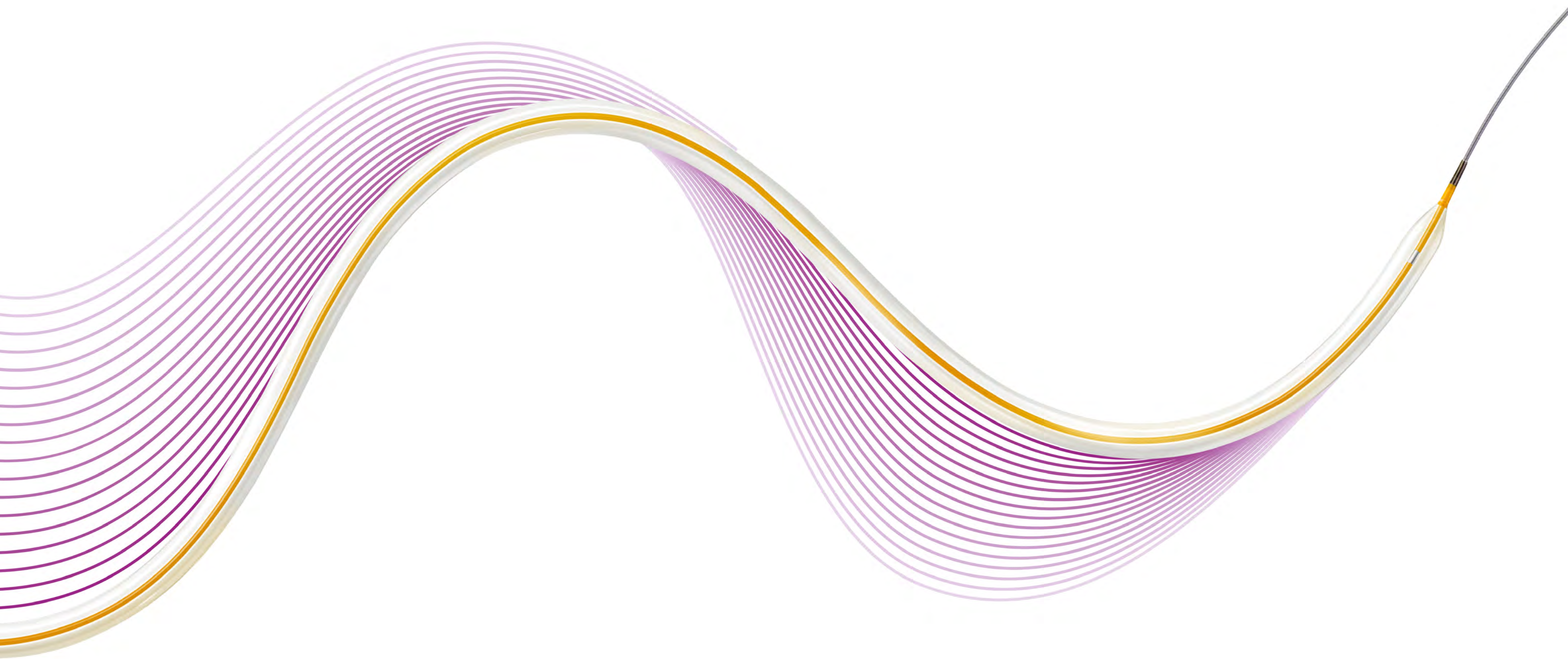


Up to 3.8 x faster deflation times¹

Due to the coaxial catheter shaft design that creates a large balloon lumen facilitating rapid inflation and deflation, Passeo-14 deflates:

3.8 x faster than Amphirion Deep

2.0 x faster than Coyote

2.4 x faster than NanoCross

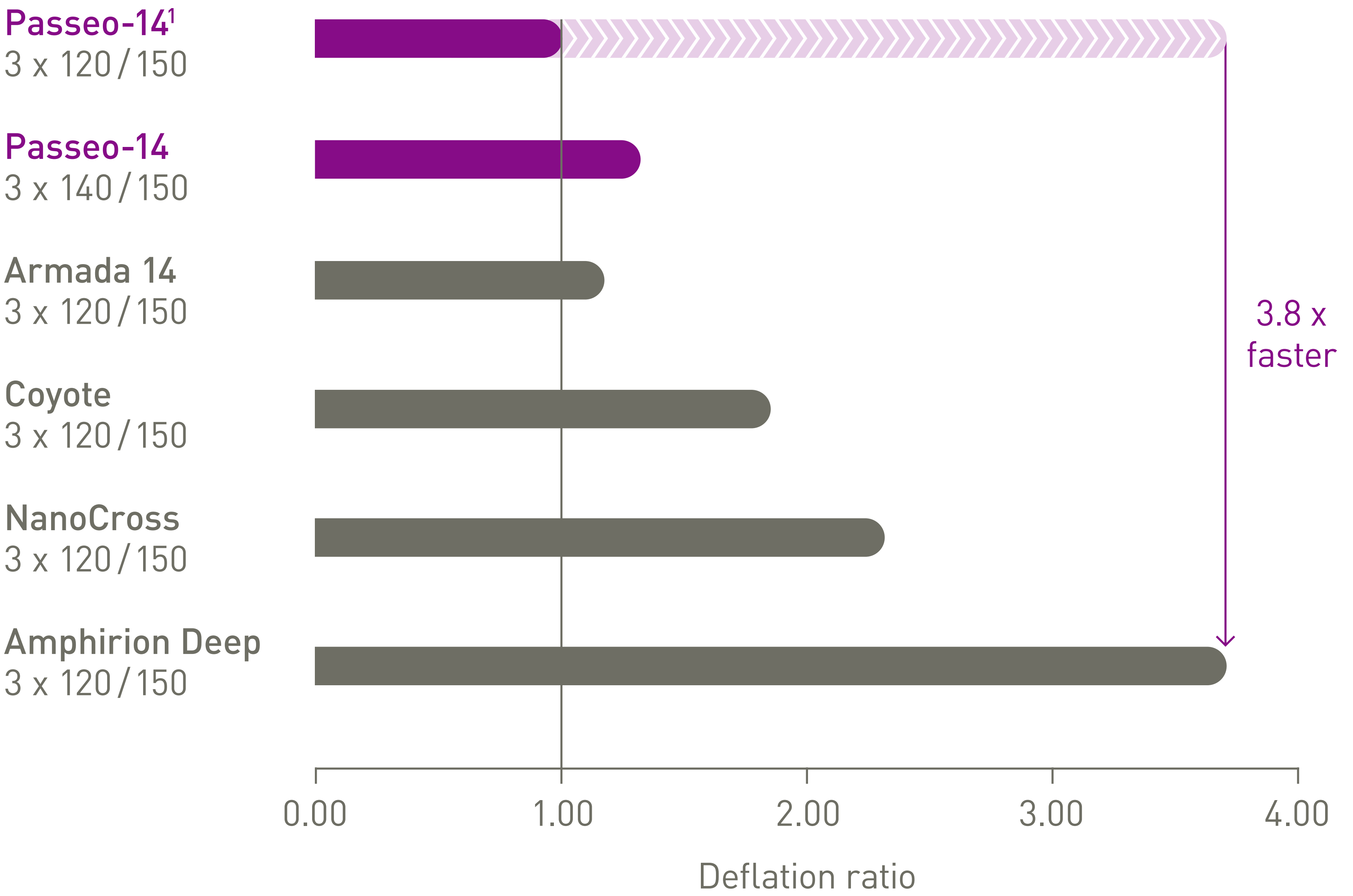

The solution for dedicated below the ankle sizes

- ø 1.5 - 2.0 mm
- 150 mm flexible distal shaft
- Tailored stiffening wire



Pre-treatment



Dilatation
pedal arch



Dilatation
plantar arch



Post dilatation


Courtesy of Dr. L. Steffanon, Vicenza, Italy






Enhanced crossability

The tri-fold balloon, which is fully coated when folded and only partly coated when inflated, enables an enhanced crossability while minimizing slippage during inflation.

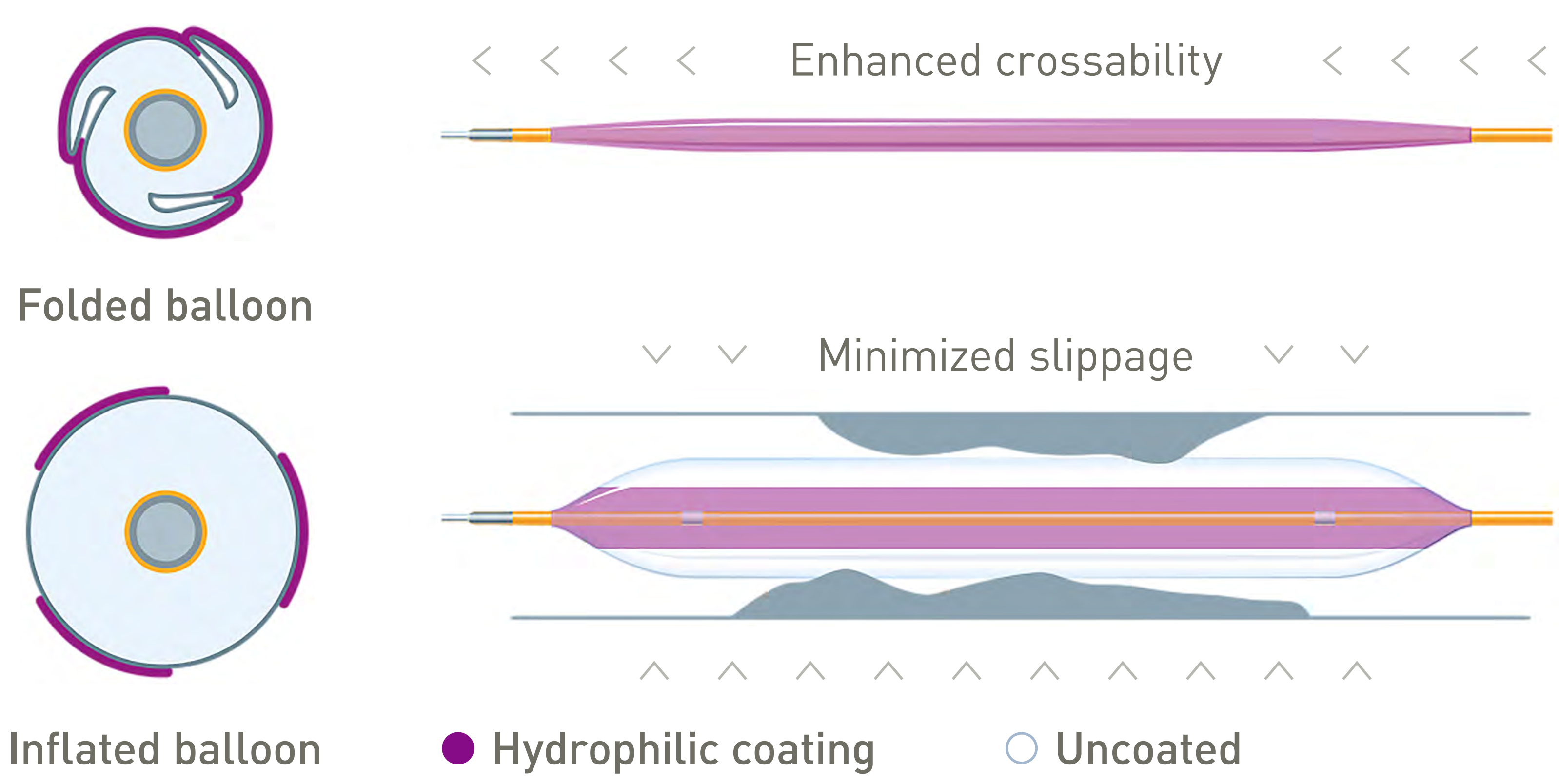


High pushability and flexibility²

Impressive pushability due to catheter shaft design featuring a unique stiffening wire in the proximal shaft of the catheter, while enabling high flexibility due to a lower profile distal shaft in small, tortuous vessels.

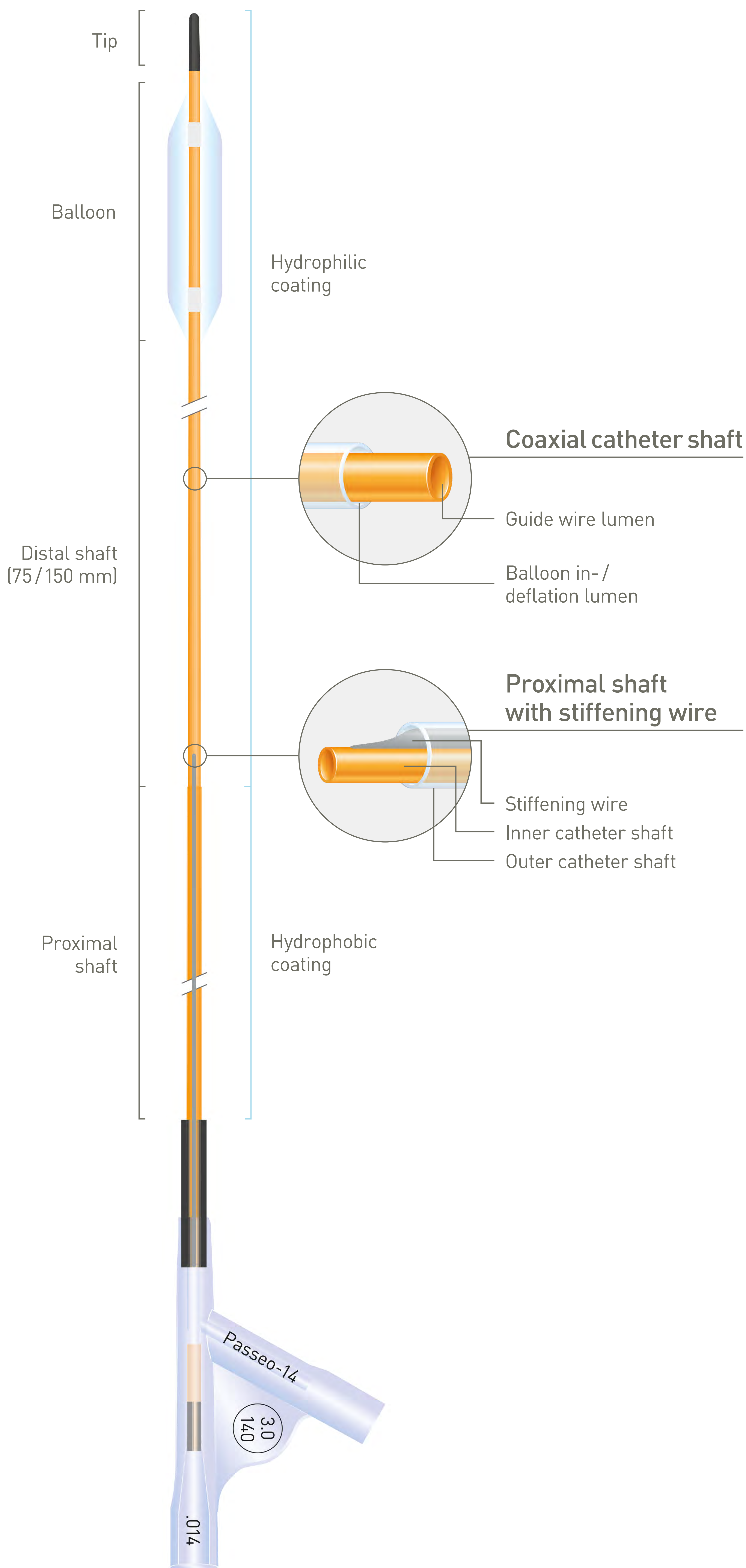

Passeo-14

Vascular
Intervention
Peripheral

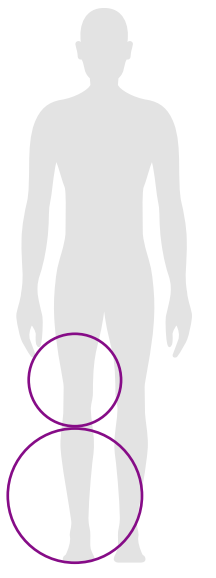

Indicated for balloon dilatation of the stenotic portion of a lower limb artery for the purpose of improving perfusion.

Technical Data	Balloon catheter
Catheter type	OTW
Recommended guide wire	0.014"
Tip	Optimized entry profile and colored
Balloon material	SCP (Semi-Crystalline Polymer), controlled compliance (4 - 6%)
Balloon folding	3-fold
Balloon coating	Hydrophilic patchwork coating
Balloon markers	2 swaged markers (zero profile)
Sizes	ø 1.5 - 4.0 mm; L: 20 - 220 mm
Distal shaft	3.1F, hydrophilic coating, coaxial design; 150 mm length (ø 1.5/2.0 x 20 - 100 mm); 75 mm length (ø 2.0 x 140 - 220 mm and ø 2.5 - 4.0 mm)
Proximal shaft	3.9F, hydrophobic coating, coaxial design; stiffening wire
Usable length	150 cm (ø 1.5 - 4.0 mm); 120 cm (ø 1.5 - 2.0 mm); 90 cm (ø 2.5 - 4.0 mm)

Compliance Chart		Balloon diameter x length (mm)					
		ø 1.5 x 20-70	ø 2.0 x 40-220	ø 2.5 x 40-220	ø 3.0 x 40-220	ø 3.5 x 40-140	ø 4.0 x 40-140
Nominal Pressure (NP)	atm*	7	7	7	7	7	7
	ø (mm)	1.5	2.0	2.5	3.0	3.5	4.0
Rated Burst Pressure (RBP)	atm*	14	14	14	14	14	14
	ø (mm)	1.57	2.08	2.61	3.18	3.63	4.16

*1 atm = 1.013 bar

Ordering Information	Catheter Length (cm)	Balloon ø (mm)	Balloon Length (mm)							
				20	40	70	100	140	180	220
4F Antegrade approach	120	1.5	380271 ^a	380277	380283	-	-	-	-	
	120	2.0	-	380278	380284	380290	380296	380302	380308	
	90	2.5	-	380279	380285	380291	380297	380303	380309	
	90	3.0	-	380280	380286	380292	380298	380304	380310	
	90	3.5	-	380281 ^a	380287 ^a	380293 ^a	380299	-	-	
	90	4.0	-	380282	380288	380294	380300	-	-	
4F Crossover approach	150	1.5	380313 ^a	380319	380325	-	-	-	-	
	150	2.0	-	380320	380326	380332	380338	380344	380350	
	150	2.5	-	380321	380327	380333	380339	380345	380351	
	150	3.0	-	380322	380328	380334	380340	380346	380352	
	150	3.5	-	380323 ^a	380329 ^a	380335 ^a	380341 ^a	-	-	
	150	4.0	-	380324	380330	380336	380342	-	-	

^a8 weeks pre-order only

1. BIOTRONIK data on file. Volume adjustment: A 3mm x 120mm balloon contains 17% less contrast media volume than a 3mm x 140mm balloon. The measured deflation time of a 3mm x 140mm balloon was adjusted by 17 % to make a direct competitive comparison; 2. BIOTRONIK Data on file.

Amphirion is a registered trademark of the Medtronic Group of Companies; Armada is a registered trademark of Abbott; Coyote is a registered trademark of Boston Scientific. NanoCross is a registered trademark of Medtronic Group of Companies.